
GB: INSTRUCTIONS FOR USE

Date of issue of the instructions for use: 01.01.2013

BENU quick pregnancy test for self-testing
Please read the instructions for use carefully before using the test!

PERFORMING THE TEST

- 1. Open the pouch shortly before testing and remove the test device.
- 2. Pull out the cap and place the absorbent tip in urine for at least 10 seconds.

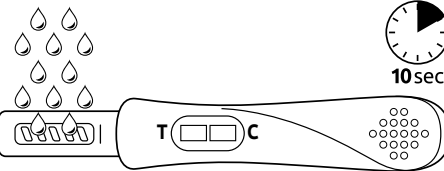Possible error: Passing the test tip with urine (e. g. in a cup) may cause malfunction of the test.

- 3. Re-cap the device and wait for 3 minutes.

READING THE RESULT

Read the result after 3 minutes. If no line appears within 3 minutes, wait 1 more minute. Do not read results after 5 minutes. A positive result may be observed in as soon as 1 minute. The colour intensity of the test bands may vary since different stages of pregnancy have different concentrations of HCG hormone.

One band:
You may assume to be **"not pregnant"**.

Two bands:
You may assume to be **"pregnant"**.

No band:
Test result **"invalid"**.

FOLLOW-UP ACTION

- If the result is **"not pregnant"** and your period late, repeat the test one day later. If the test is still negative, consult your gynaecologist.
- If the result is **"pregnant"**, please consult your gynaecologist to confirm your pregnancy (avoid potential risks to harm the embryo such as e. g. x-ray).
- If the result is **"invalid"**, the absorbent tip probably was not thoroughly wet. Please repeat the testing with a new test.

If you repeat the test several times and results are contradictory, please consult your gynaecologist.

QUESTIONS & ANSWERS

- 1. **How does the BENU QUICK pregnancy test work?** Human Chorionic Gonadotropin (hCG) is a hormone, produced by the developing placenta shortly after the conception and secreted into the urine. It increases up to the 3rd month of pregnancy. The BENU QUICK pregnancy test detects this hormone in urine rapidly.
- 2. **When to do the test?** Any time of day from missed period.
- 3. **Can alcohol alter the test result?** No, alcohol will usually not have any effect on the test result.
- 4. **Can medication alter the test result?** hCG, FSH, LH and TSH based medication may have an effect on the test result, other hormone-based contraceptives like the birth control pill may not.

LIMITATIONS

- Test result must be read in 5 minutes. A negative result, which may change to positive after 5 minutes, is not an accurate reading.
- Due to the high sensitivity of the test, even very low hCG values can be verified in the early stages of pregnancy. During this early stage, the number of spontaneous abortions is about 31 %, approx. 22 % remain clinically undetected.
- Beside pregnancy, some illnesses and spontaneous abortions may lead to increased levels of hCG in urine, which may be detected by the test.
- Menopausal urine may similarly cause false positives.
- Low concentration of hCG in a very early stage of pregnancy can give a negative result. In this case, another specimen should be obtained at least 48 hours later and tested.
- Sterility treatments, based on hCG, may cause false results. Please consult your gynaecologist.

PRECAUTIONS

- Keep out of the reach of children, do not swallow!
- Do not disassemble the test device!
- Discard the test immediately after use!
- Do not use the test if the pouch is damaged or beyond expiry date!

PRODUCT INFORMATION

PRODUCT CODE:	BN-513
PACKAGE CONTENT:	1 test, 1 instructions for use.
PRODUCT CODE:	BN-520
PACKAGE CONTENT:	2 tests, 1 instructions for use.
STORAGE:	At room temperature (4°C to 30°C) away from direct sunlight, moisture and heat.
CHARACTERISTICS:	Sensitivity from 25 IU/L. Accuracy over 99%.
ACTIVE INGREDIENTS:	Monoclonal antibodies against hCG.

MEANING OF SYMBOLS

- Follow the instructions for use
- Tests per package
- Store between 4°C and 30°C
- Manufacturer
- In vitro use. Not to be swallowed!
- Expiry date
- Single use only
- Batch code

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